

Hand-bikes by

COSTA hand-bike is a light and versatile 3-wheel bike, equipped with a 27 or 8 gear shift suitable for all terrains. A special aid for running on the free time. The high seat and the 90° leg position make the bike unique for speed and easiness of transfers and the wide range of adjustments satisfies every user needs.

CE



COSTA



Seat height adjustment

Foot-rest depth adjustment



The seat tube can be fixed under or above the wheel axle to increase or to reduce the rear seat/floor height.

Hand-bikes

COSTA

Technical features

- Aluminium frame.
- 26" wheels assembled with pneumatic tyres.
- Height, depth, tilt adjustable seat and back-rest.
- Foot-rest adjustable in several positions.
- Rear wheel cambers: 12° or 6°.
- Brake and shift controls on handles.
- Twin front brake.
- Parking brake.
- Tension adjustable back-rest upholstery.
- Leg keeper and foot-holding belts.
- Rear detachable bumper bar.
- Wireless computer.
- Weight: 16 kg (standard set-up).



Shift handle



Disk brake



Bumper bar



Front light

COSTA 8CR

Revo 8 gear hub shift and coaster brake.

COSTA 8F

8 gear hub shift with lever brake.

COSTA 27

Shimano XT 27 gear shift group with 22-32-44 chainwheel teeth cranksset, supplied with a 9 cassette sprockets:
11-12-14-16-18-21-24-28-32 or
12-13-14-15-16-17-19-21-23 sets.



F




H

Tetraplegic handles

Order form

COSTA 27 COSTA 8CR COSTA 8F

- Seat width: 36 38 40 42 cm
- Seat depth: 38 40 cm
- Back-rest height: cm
- Distance crankset centre/back-rest: (A) cm
- Distance seat upholstery/crankset centre: (B) cm
- Front seat height from the foot-rest: (C) cm (from 33 to 42 cm).
- Rear seat height from the floor: 34,5 cm 38,5 cm
- Frame colour: Ferrari red metallic grey
- Cassette sprockets (27 gear only): 11-12-14-16-18-21-24-28-32
 12-13-14-15-16-17-19-21-23
- Rear wheel camber: 6° 12°

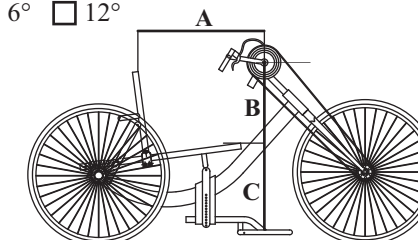

Standard measurements:

- A = 55 cm
- B = 40 cm
- C = 40 cm

Options

- Spoke guards: with **OFFCARR** logo without **OFFCARR** logo
- Bell: Rear-view mirror:
- Front light: Disk brake:
- Lined foam cushion: 3 cm 5 cm 7 cm
- Water-bottle holder: Water-bottle: Termic water-bottle:
- Padded rigid back-rest:
- Off-road wheels: Flag: Bumper bar:
- Other frame colour: (see "Frame colour" leaflet)
- Tetraplegic handles: model (F) model (H)



Via dell' Artigianato II, 29
35010 Villa del Conte (PD) Italy
Tel. +39 049 9325733 - Fax +39 049 9325734
www.offcarr.com - offcarr@offcarr.com

Dealer address

Date: _____
Order number: _____
Custume ref.: _____
Notes: _____